


# Piccolo<sup>®</sup> Lipid Panel

**For In Vitro Diagnostic Use and Healthcare Professional Use Only**  
**Customer and Technical Service: 1- 800-822-2947**  
**Customers outside the US: +49 6155 780 210**  
[AB-PiccoloTechSupport@zoetis.com](mailto:AB-PiccoloTechSupport@zoetis.com)

**Applicable to US customers only**  
**CLIA Waived: Use lithium heparin whole blood, only**  
**Moderate Complexity: Use lithium heparin whole blood, lithium heparin plasma, or serum**



Abaxis, Inc.  
 3240 Whipple Rd.  
 Union City, CA 94587  
 USA

CH	REP
----	-----

Tao of Excellence GmbH  
 Vorstadt 26,  
 8200 Schaffhausen  
 Switzerland

EC	REP
----	-----

ABAXIS Europe GmbH  
 Bunsenstr. 9-11  
 64347 Griesheim  
 Germany

## 1. Intended Use

The Piccolo<sup>®</sup> Lipid Panel reagent disc, used with the Piccolo Xpress<sup>®</sup> Chemistry analyzer, is intended for the *in vitro* quantitative determination of total cholesterol (CHOL), high-density lipoprotein cholesterol (HDL), and triglycerides (TRIG) in capillary (fingerstick) lithium heparinized whole blood, venous lithium heparinized whole blood, lithium heparinized plasma, or serum in a clinical laboratory setting or point-of-care location/ near- patient testing. From these determinations low-density lipoprotein cholesterol (LDL), very low-density lipoprotein cholesterol (VLDL), non-HDL cholesterol, and a total cholesterol/high-density lipoprotein cholesterol ratio (TC/H) are calculated by the analyzer.

Lipid measurements are used in the diagnosis and treatment of lipid and lipoprotein disorders, atherosclerosis, various liver and renal diseases, diabetes mellitus, and other diseases involving lipid metabolism or various endocrine disorders.

### For US Customers Only

This test is waived under CLIA '88 regulations. If a laboratory modifies the test system instructions, then the test is considered high complexity and subject to all CLIA requirements. For CLIA waived labs, only lithium heparin whole blood (venous or capillary) may be tested. For use in moderate complexity labs, lithium heparinized whole blood, lithium heparinized plasma, or serum may be used.

A CLIA Certificate of Waiver is needed to perform CLIA waived testing. A Certificate of Waiver can be obtained from the Centers for Medicare & Medicaid Services (CMS).

## 2. Summary and Explanation of Tests

### Clinical Significance

In adults, measurement of the serum lipids and lipoproteins is useful in characterizing an individual's risk of developing cardiovascular diseases (CVD).<sup>1</sup> Consensus-based guidelines for measurement and cut-points for interpretation have been provided by the National Cholesterol Education Program.<sup>2,3,4</sup>

The circulating lipids are carried on lipoproteins. The LDL fraction, the major lipoprotein contributor to the development of atherosclerosis and for which treatment has been conclusively demonstrated to be effective, carries most of the circulating cholesterol in the blood. Total serum cholesterol has been measured for many years to quantify the total amount of lipoproteins as a convenient means of assessing CVD risk. However, some of the cholesterol is carried on HDL particles, which are anti-atherogenic or inversely associated with risk of developing CVDs. Thus, quantitation of the major individual cholesterol-carrying lipoproteins, LDL and HDL, provides a better assessment of overall risk.

Triglycerides, the body's major fuel, are carried into the blood stream on large lipoproteins called chylomicrons (CM). VLDL particles also carry triglycerides, primarily synthesized in the liver from excess fatty acids. In the circulation triglycerides are hydrolyzed and their fatty acids transported into peripheral cells leaving remnant particles, precursors to LDL. After an overnight fast, chylomicrons have generally been cleared from the circulation. Higher levels of triglycerides measured in a fasting specimen indicate impaired clearance or over-production, which may increase risk of developing CVD, making their measurement useful in characterizing metabolic disorders and overall risk.

The US National Heart, Lung and Blood Institute organized the National Cholesterol Education Program, which convened

expert panels to develop clinical guidelines for classification and treatment of high cholesterol. The most recent recommendations, the Adult Treatment Panel III guidelines,<sup>2,3,4</sup> base treatment decisions primarily on the LDL levels, calculated as part of the lipid panel after measurement of total cholesterol, HDL, and triglycerides. LDL cut-points of 100, 130, 160, and 190 mg/dL (2.59, 3.37, 4.14 and 4.92 mmol/L) define optimal, near optimal, borderline high, high and very high risk categories. An HDL value below 40 mg/dL (1.04 mmol/L) is low, considered to be a risk factor by the ATP III, modifying the LDL treatment goal. An HDL value above 60 mg/dL (1.55 mmol/L) is defined as high, considered desirable and a negative risk factor, subtracting from the total number of risk factors in selecting the appropriate treatment goal for LDL. For triglycerides, cut-points of 150, 200, and 500 mg/dL (1.7, 2.26 and 5.65 mmol/L) define normal, borderline-high, high, and very-high levels. Additionally, the calculated non-HDL cholesterol optimally should be < 130 mg/dL (3.37 mmol/L), with increased risk for cardiovascular disease associated with concentrations of 130 – 189 mg/dL (3.37 – 4.9 mmol/L), and high risk of CVD associated with values > 189 mg/dL (4.9 mmol/L).


The Piccolo Total Cholesterol and HDL assays are not certified by the Cholesterol Reference Method Laboratory Network (CRMLN), which is coordinated by the Centers for Disease Control. The requirements of accuracy and precision have been previously established to meet the accuracy and precision requirements of CRMLN.

**As with any diagnostic test procedure, all other test procedures including the clinical status of the patient, should be considered prior to final diagnosis.**

### 3. Principles of Methods

#### Total Cholesterol (CHOL)





The Abaxis CHOL assay is an enzymatic end-point method that uses cholesterol esterase (CE) and cholesterol dehydrogenase (CHDH).<sup>5</sup>

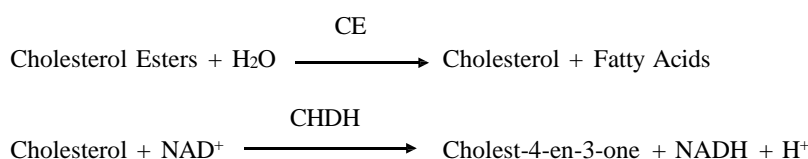

CE hydrolyzes cholesterol esters to form cholesterol and fatty acids. The CHDH reaction converts cholesterol to cholest-4-en-3-one. The NADH is measured bichromatically at 340 nm and 405 nm. NADH production is directly proportional to the amount of cholesterol present. An assay-specific blank is also monitored to ensure no extraneous reactions interfere with the calculations of CHOL levels.

#### High-Density Lipoprotein Cholesterol (HDL)


The Abaxis HDL assay is a precipitation method that utilizes polyethylene glycol-modified cholesterol esterase (PEG-CE) and cholesterol oxidase (PEG-CO) for additional specificity.<sup>6</sup> The reaction mechanism follows:

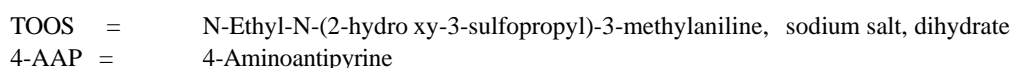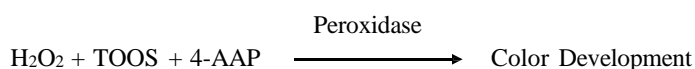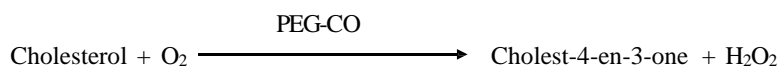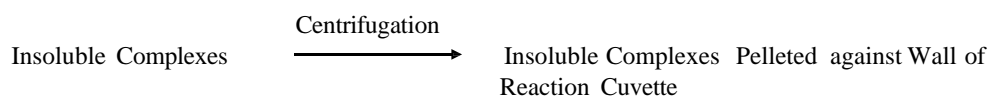

The precipitating agents dextran sulfate and magnesium sulfate (MgSO<sub>4</sub>) specifically form insoluble complexes with chylomicrons (CM), VLDL, and LDL in plasma or serum. The insoluble complexes are pelleted to the wall of the reaction cuvette within the analyzer. The remaining HDL is hydrolyzed by PEG-CE to make cholesterol and fatty acids. Cholesterol reacts with PEG-CO to produce cholest-4-en-3-one and peroxide (H<sub>2</sub>O<sub>2</sub>). The peroxidase reaction results in the production of a purple colored product that has an absorbance maximum at 550 nm and is referenced to absorbance at 630 nm. HDL cholesterol concentration is directly proportional to the absorbance maximum in this end-point reaction. A sample blank is also monitored to ensure no extraneous reactions interfere with the calculations of HDL levels.

### Triglycerides (TRIG)

The Abaxis TRIG assay is an enzymatic end-point method that makes use of four enzymes.<sup>7,8</sup> The reaction mechanism follows:

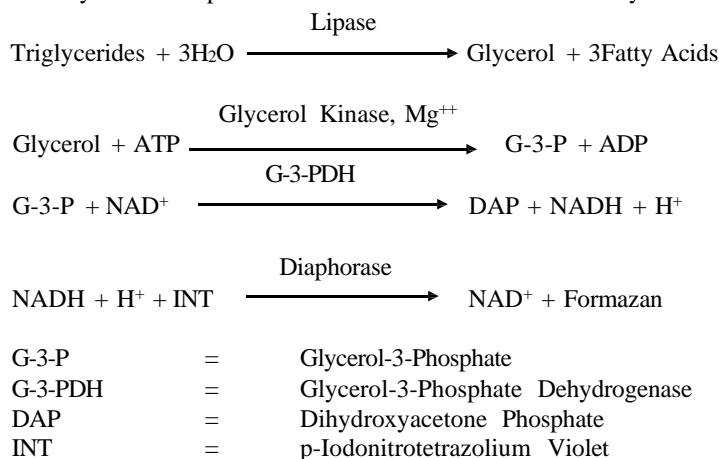

In the first step, the triglycerides are hydrolyzed into glycerol and fatty acids in a reaction catalyzed by lipase. Glycerol is then phosphorylated in an ATP-requiring reaction catalyzed by glycerol kinase (GK). The glycerolphosphate is then oxidized to dihydroxyacetone phosphate with the simultaneous reduction of NAD<sup>+</sup> to NADH in a reaction catalyzed by glycerol-3-phosphate dehydrogenase (G-3-PDH). The NADH is then oxidized with the simultaneous reduction of INT in a reaction catalyzed by diaphorase. The intensity of the highly colored formazan is measured bichromatically at 500 nm and 850 nm and is directly proportional to the concentration of triglycerides in the sample. An assay-specific blank is also monitored to ensure no extraneous reactions interfere with the calculations of TRIG levels. The results provide a measure of total triglycerides without a glycerol blank.

### LDL (Calculated)

The Piccolo analyzer or Piccolo Xpress analyzer automatically calculates the concentration of LDL in each sample using the directly determined values for total cholesterol, HDL, and triglycerides and the standard Friedewald equation.<sup>9</sup> This equation is not valid for triglyceride concentrations above 400 mg/dL (4.52 mmol/L), non-fasting patients, and patients with Type III hyperlipoproteinemia (dysbetalipoproteinemia).<sup>9,10</sup> An LDL value is not reported for samples with triglycerides greater than 400 mg/dL (4.52 mmol/L) or if any of the directly measured analyte values is unavailable. On the printout, the calculated value for LDL is followed by a “C” to indicate that it is calculated.

### VLDL (Calculated)

The Piccolo analyzer or the Piccolo Xpress analyzer automatically calculates the concentration of VLDL in each sample using the standard triglycerides/5 (if units in mg/dL) equation.<sup>9</sup> This equation is not valid for triglyceride concentrations above 400 mg/dL (4.52 mmol/L), non-fasting patients, and patients with Type III hyperlipoproteinemia (dysbetalipoproteinemia).<sup>9,10</sup> Of course, no VLDL value is calculated if no triglyceride value is available. On the printout, the calculated value for VLDL is followed by a “C” to indicate that it is calculated.

### Total Cholesterol/HDL Ratio (Calculated)

The Piccolo analyzer or the Piccolo Xpress analyzer automatically calculates the total cholesterol/HDL ratio (abbreviated as TC/H) for each sample. If the directly measured total cholesterol or HDL value is missing, no ratio is provided. On the printout, the calculated value for TC/H is followed by a “C” to indicate that it is calculated.

### Non-HDL (Calculated) – only available on the Xpress

The Piccolo Xpress analyzer automatically calculates the non-HDL cholesterol (abbreviated as nHDLc) for each sample. The nHDLc is calculated by subtracting HDL cholesterol from total cholesterol (CHOL). On the printout, the calculated value for nHDLc is followed by a “C” to indicate that it is calculated. If the directly measured total cholesterol (CHOL) or HDL value is missing, nHDLc is not calculated.

## 4. Principle of Operation

Refer to the Piccolo Xpress chemistry analyzer Operator’s Manual, for the Principles and Limitations of the Procedure.

## 5. Description of Reagents

### Reagents

Each Piccolo Lipid Panel reagent disc contains dry test-specific reagent beads (described below). A dry sample blank reagent bead (comprised of buffer, surfactants, excipients, and preservatives) is included in each disc for use in calculating concentrations of HDL. Dedicated blank beads are also included in the disc to calculate concentrations of CHOL and TRIG. Each disc also contains a diluent consisting of a surfactant and preservatives.

**Table 1: Reagents**

Component	Quantity/Disc
4-Aminoantipyrine	12 µg
Adenosine 5'-triphosphate, disodium salt	9.2 µg
Ascorbate oxidase	0.042 U
Cholesterol dehydrogenase	0.27 U
Cholesterol esterase	0.27 U
Dextran sulfate	8.4 µg
Diaphorase	0.25 U
3,5-Dichloro-2-hydroxybenzenesulfonic acid, Sodium salt (DHBSA)	34.8 µg
Ethylenediaminetetraacetic acid (EDTA), Tetrasodium Salt	180 µg
N-Ethyl-N-(2-hydroxy-3-sulfopropyl)-3-methylaniline, sodium salt, dihydrate (TOOS)	79 µg
Glycerol kinase	0.084 U
Glycerol-3-phosphate dehydrogenase	0.16 U
Iodonitrotetrazolium chloride (INT)	8.4 µg
Lipase	16.8 U
Magnesium chloride, hexahydrate	8.6 µg
Magnesium sulfate, heptahydrate	197 µg
Nicotinamide adenine dinucleotide, monosodium salt (NAD)	455 µg
PEG-cholesterol esterase	0.0034 U
PEG-cholesterol oxidase	0.1280 U
Peroxidase	0.03 U
Uricase	0 U
Buffers, surfactants, excipients, and preservatives	

### Warnings and Precautions

- For *In vitro* Diagnostic Use
- The diluent container in the reagent disc is automatically opened when the analyzer drawer closes. A disc with an opened diluent container cannot be re-used. Ensure that the sample or control has been placed into the disc before closing the drawer.
- Used reagent discs contain human body fluids. Follow good laboratory safety practices when handling and disposing of used discs.<sup>11</sup> See the Piccolo Xpress chemistry analyzer Operator’s Manual for instructions on cleaning biohazardous spills.

- The reagent discs are plastic and may crack or chip if dropped. Never use a dropped disc as it may spray biohazardous material throughout the interior of the analyzer.
- Reagent beads may contain acids or caustic substances. The operator does not come into contact with the reagent beads when following the recommended procedures. In the event that the beads are handled (e.g., cleaning up after dropping and cracking a reagent disc), avoid ingestion, skin contact, or inhalation of the reagent beads.
- This mixture does not contain substances assessed to be vPvB / PBT according to Regulation (EC) No 1907/2006, Annex XIII. This product contains components considered to have endocrine disrupting properties for environment, according to REACH Article 57(f), Regulation (EU) 2018/605 or Regulation (EU) 2017/2100. Refer to the Safety Data Sheet (SDS) for specific details.
- Materials incorporated into the device may consist of substances derived from other microbial and animal sourced materials.

### **Instructions for Reagent Handling**

Reagent discs may be used directly from the refrigerator without warming. A disc not used within 20 minutes of opening the pouch should be discarded. Do not allow discs sealed in their foil pouches to remain at room temperature at 20-25°C (68-77°F) longer than 48 hours prior to use. Open the sealed foil pouch, remove the disc, and use according to the instructions provided in the Piccolo Xpress chemistry analyzer Operator's Manual.

### **Storage**

Store reagent discs in their sealed pouches at 2-8 °C (36-46 °F). Do not expose opened or unopened discs to direct sunlight or temperatures above 32 °C (90 °F). Reagent discs may be used until the expiration date printed on the box label and on each pouch. The expiration date is also encoded in the bar code printed on the bar code ring. An error message will appear on the Piccolo Xpress chemistry analyzer display if the reagents have expired.

### **Indications of Reagent Disc Instability/Deterioration**

A torn or otherwise damaged pouch may allow moisture to reach the unused disc and adversely affect reagent performance. Do not use a disc from a damaged pouch.

## **6. Instrument**

See the Piccolo Xpress chemistry analyzer Operator's Manual for complete information on use of the analyzer.

## **7. Sample Collection and Preparation**

Sample collection techniques are described in the "Sample Collection" section of the Piccolo Xpress chemistry analyzer Operator's Manual.

- According to the ATP III,<sup>2,3,4</sup> fasting samples (eight to 12 hours) should be used to determine CHOL, HDL, TRIG, and LDL. Hence, it is highly recommended that fasting samples be used with the Lipid Panel Disc. Should the patient not be fasting, the TRIG and calculated LDL values are not valid.
- The minimum required sample size is ~100 µL for all sample types and controls. The reagent disc sample chamber can contain up to 120 µL of sample.
- Whole blood venipuncture samples should be run within 60 minutes of collection.<sup>12,13</sup> The sample may be separated into plasma or serum and stored in capped sample tubes at 2-8 °C (36-46 °F) if the sample cannot be run within 60 minutes.
- Start the test within 10 minutes of transferring the sample into the reagent disc.

### **Venous Sampling Notes:**

- Use only lithium heparin evacuated specimen collection tubes with or without gel separators for whole blood or plasma samples. Use no-additive evacuated specimen collection tubes or serum separator tubes for serum samples.
- Whole blood samples obtained by venipuncture must be homogeneous before transferring a sample to the reagent disc. Gently invert the collection tube several times just prior to sample transfer. Do not shake the collection tube; shaking may

cause hemolysis.

### **Capillary Sampling Notes:**

- Proper and adequate capillary sampling techniques are required for reliable results.
- For capillary blood collection, the patient's hands should be washed thoroughly with soap (without glycerin or glycerol) and warm water. Assure complete drying.
- For capillary sampling, gentle squeezing along the finger up to the puncture site is acceptable, but do not milk the finger to get the blood to flow, which can lead to excess hemolysis.

## **8. Procedure**

### **Materials Provided**

- One Piccolo Lipid Panel PN: 400-1025 (a box of discs PN: 400-0025)

### **Materials Required but not Provided**

- Piccolo Xpress chemistry analyzer
- A sample transfer pipette (fixed volume approximately 100 µL) and tips are provided with each Piccolo Xpress chemistry analyzer and may be reordered from Abaxis.
- Commercially available control reagents recommended by Abaxis (contact Abaxis Technical Service for approved control materials and expected values).
- Timer

### **Test Parameters**

The Piccolo Xpress chemistry analyzer operates at ambient temperatures between 15°C and 32 °C (59-90 °F). The analysis time for each Piccolo Lipid Panel Reagent Disc is less than 15 minutes. The analyzer maintains the reagent disc at a temperature of 37 °C (98.6 °F) over the measurement interval.

### **Test Procedure**

The complete sample collection and step-by-step operating procedures are detailed in the Piccolo Xpress chemistry analyzer Operator's Manual.

### **Calibration**

The Piccolo Xpress chemistry analyzer is calibrated by the manufacturer before shipment. The bar code printed on the bar code ring provides the analyzer with disc-specific calibration data. See the Piccolo Xpress chemistry analyzer Operator's Manual.

The Piccolo Xpress chemistry analyzer is calibrated using internal calibrators or reference materials. The methods and materials used in the procedures for the calibration and control value assignments are traceable either to the comparative method and/or to the standards listed below:

- IFCC – International Federation of Clinical Chemistry
- NIST – National Institute of Standards and Technology
- CRMLN – Cholesterol Reference Method Laboratory Network
- SRM – Standard Reference Material
- CRM – Certified Reference Material

The calibration and value assignment processes are in compliance with EN ISO 17511, Metrological Traceability of Values Assigned to Calibrators and Control Materials.

The measured analytes in the Lipid Panel are traceable to the following reference materials or methods. The Piccolo Xpress Analysis System controls and calibration verification materials are validated for use only with the Piccolo Xpress Analysis System and assigned values may not be commutable with other methods. For specific analyte reference standard methods, refer to Section 3. Principles of Method.

Analyte	Calibration Standard	Method
CHOL	CRMLN	Enzymatic
HDL	CRMLN	Enzymatic
TRIG	Roche Cobas Fara	Enzymatic
LDL	N/A	Calculated from CHOL, HDL, and TRIG results using the Friedewald Equation.

### Quality Control

See Section 6 (Calibration and Quality Control) of the Piccolo Xpress Operator's Manual. Performance of the Piccolo Xpress chemistry analyzer can be verified by running controls. For a list of approved quality control materials with acceptance ranges, please contact Abaxis Technical Support. Other human serum or plasma-based controls may not be compatible. Quality control materials should be stored as per the package -insert included with the controls.

If control results are out of range, repeat one time. If still out of range, call Abaxis Technical Support. Do not report results if controls are outside their labeled limits. See the Piccolo or Piccolo Xpress Operator's Manual for a detailed discussion on running, recording, interpreting, and plotting control results.

**Waived Laboratories (US only):** Abaxis recommends control testing as follows:

- at least every 30 days
- whenever the laboratory conditions have changed significantly, e.g. Piccolo moved to a new location or changes in temperature control
- when training or retraining of personnel is indicated
- with each new lot (CLIA waived tests in waived status labs)

**Non-Waived Laboratories (US only):** Abaxis recommends control testing to follow federal, state, and local guidelines.

And the non-waived ones will use the below:

- Performance of the Piccolo Xpress chemistry analyzer can be verified by running controls. For a list of approved quality control materials with acceptance ranges, please contact Abaxis Technical Support. Other human serum or plasma-based controls may not be compatible. Quality control materials should be stored as per the package-insert included with the controls.
- If control results are out of range, repeat one time. If still out of range, call Abaxis Technical Support. Do not report results if controls are outside their labeled limits. See the Piccolo Xpress chemistry analyzer Operator's Manual for a detailed discussion on running, recording, interpreting, and plotting control results.

## 9. Results

The Piccolo Xpress chemistry analyzer automatically tests, calculates, and prints the analyte concentrations in the sample. For those results that are calculated and not directly determined, LDL, VLDL, and TC/H, the results are followed by a "C" to indicate that they are calculated. Details of the endpoint and rate reaction calculations are found in the Piccolo Xpress chemistry analyzer Operator's Manual.

Interpretation of results is detailed in the Operator's Manual. Results are printed onto result cards or paper rolls available from Abaxis. The result cards or paper rolls have an adhesive backing for easy placement in the patient's files.

The reaction for each analyte occurs at 37°C (98.6°F).

## 10. Limitations of Procedure

General procedural limitations are discussed in the Piccolo Xpress chemistry analyzer Operator's Manual.

- The only anticoagulant **recommended for use** with the Piccolo Xpress chemistry analyzer is **lithium heparin**. Abaxis has performed studies demonstrating that EDTA, fluoride, oxalate, and any anticoagulant containing ammonium ions will interfere with at least one chemistry contained in the Piccolo Lipid Panel. Do not use sodium heparin.
- Samples with hematocrits in excess of 62% packed red cell volume may give inaccurate results. Samples with high hematocrits may be reported as hemolyzed. These samples may be spun down to get plasma and then re-run in a new reagent disc.
- **Any result for a particular test that exceeds the assay range should be analyzed by another approved test method or sent to a referral laboratory. Do not dilute the sample and run it again on the Piccolo Xpress chemistry analyzer.**

**Warning:** Extensive testing of the Piccolo Xpress chemistry system has shown that, in very rare instances, sample dispensed into the reagent disc may not flow smoothly into the sample chamber. Due to the uneven flow, an inadequate quantity of sample may be analyzed and several results may fall outside the expected ranges. The sample may be re-run using a new reagent disc.

**Note:** Operators shall report any serious incident that has occurred in relation to the device to the manufacturer.

### Interference

Substances were tested as interferents with the analytes. Human serum pools were prepared. The concentration at which each potential interferent was tested was based on the testing levels in CLSI EP7-A.<sup>15</sup>

### Effects of Endogenous Substances

- Physiological interferents (hemolysis, icterus and lipemia) cause changes in the reported concentrations of some analytes. The sample indices are printed on the bottom of each printout to inform the operator about the levels of interferents present in each sample.
- The Piccolo Xpress chemistry system suppresses any results that are affected by >10% interference from hemolysis, lipemia or icterus. “HEM”, “LIP”, or “ICT” respectively, is printed on the printout in place of the result.
- For maximum levels of endogenous substances contact Abaxis Technical Service.

### Effects of Exogenous and Therapeutic Substances

Sixteen therapeutic substances were selected as potential interferents for the total cholesterol, HDL, and triglycerides assays based on recommendations by Young.<sup>16</sup> Significant interference is defined as a >10% shift in the result for a normal range specimen. Human serum pools were supplemented with a known concentration of the drugs or chemicals and then analyzed.

**Please see Table 2 for a list of exogenous and therapeutic substances evaluated. Please see Table 3 for a list of analytes where interference was observed.**

**Table 2: Therapeutic Substances Evaluated**

	Highest Concentration Tested (mg/dL)
Acetaminophen	100
Acetoacetate, Lithium	102
Acetylsalicylic Acid	50
Ascorbic Acid	3
Digoxin	5
Glutathione	30
Heparin, Lithium	4.4 (800 U/dL)
β-Hydroxybutyrate	1,000
Ibuprofen	50
Isoniazide	4
Lactate, Lithium	84
Lidocaine	1
Methicillin, Sodium	100
Phenytoin	3



Salicylic Acid  
Theophylline

50  
20

None of the assays demonstrated significant interference at the concentrations of exogenous substances tested.

**Table 3: Concentration of Analytes in Serum Pools (2 Levels) Used for Exogenous Substance Interference Studies**

Analyte	Concentration
<b>Cholesterol (CHOL)</b>	163 and 197 mg/dL (4.22 and 5.1 mmol/L)
<b>HDL</b>	39 and 52 mg/dL (1.01 and 1.35 mmol/L)
<b>Triglycerides (TRIG)</b>	125 and 183 mg/dL (1.41 and 2.08 mmol/L)

## 11. Cutpoints

Consensus based cutpoints for the lipid/lipoprotein analytes have been established by the National Cholesterol Education Program as follows:<sup>2,3,4</sup>

**Table 4: Medical Decision Values<sup>2,3,4</sup>**

	Interpretation	Cutpoints mg/dL	Cutpoints mmol/L
Total Cholesterol (CHOL)	Desirable	< 200	< 5.17
	Borderline High	200–239	5.17–6.18
	High	≥ 240	6.20
HDL	Low HDL - Risk Factor	< 40	< 1.03
	High HDL - Negative Risk Factor (Desirable)	≥ 60	≥ 1.55
Triglycerides (TRIG)	Normal	< 150	<1.70
	Borderline High	150–199	1.70–2.25
	High	200–499	2.26–5.64
	Very High	≥ 500	≥ 5.65
LDL*	Optimal	< 100	< 2.58
	Near Optimal	100–129	2.58–3.33
	Borderline High	130–159	3.36–4.11
	High	160–189	4.13–4.88
	Very High	≥ 190	≥ 4.91
VLDL (CALC) **	Normal	< 30	
	High	≥ 30	
nHDLc (CALC)**	Optimal	< 130	< 3.37
	Increased Risk	130–189	3.37–4.90
	High Risk	> 189	> 4.90
Total Chol/HDL Ratio	Low Risk	≤ 5	≤ 4.5
	High Risk	> 5	> 4.5
		Male	Female

\* The Piccolo or Piccolo Xpress calculates the LDL concentration using the Friedewald equation.<sup>9</sup>

\*\* For further information see: NCEP, ATP III Report 2002, Section II. Rationale for Intervention, 3. Other Lipid

**Total Cholesterol / HDL Ratios (TC/H)**

The total cholesterol to HDL ratio (TC/H) is calculated as a convenience to users. A TC/H ratio of  $\leq 5$  is generally considered desirable for men. Because women usually have higher HDL values than men, some recommend a cutpoint of 4.5 for women.<sup>14</sup> This ratio has been advocated by some as a simple and convenient means of expressing CVD risk in a single number, incorporating total cholesterol as a marker for atherogenic lipoproteins in the numerator and the anti-atherogenic HDL cholesterol in the denominator.<sup>1</sup> User's should be aware that even though the TC/H ratio is a powerful predictor of CVD risk as shown by numerous epidemiology studies,<sup>1</sup> the NCEP does not recommend its use in managing patients. The NCEP clinical guidelines base treatment decisions on the individual lipoproteins (Table 5) and consider use of the ratio as a possible diversion from the priority, the individual lipoprotein measurements.<sup>2,3,4</sup>

**Non-HDL Cholesterol (nHDLc)**

NCEP, ATP III Report 2002, reported the clinical utility of the calculated nHDLc. Studies have demonstrated that non-HDL cholesterol shows a stronger correlation with coronary mortality when compared with LDL cholesterol. Moreover, non-HDL cholesterol is highly correlated with total apolipoprotein B (apoB), the primary apolipoprotein associated with all atherogenic lipoproteins.”<sup>2</sup>

**12. Performance Characteristics**

**Linearity**

The chemistry for each analyte is linear over the dynamic range listed below when the Piccolo Xpress chemistry analyzer or the Piccolo Xpress chemistry analyzer is operated according to the recommended procedure (refer to Piccolo Xpress chemistry analyzer Operator’s Manual). This evaluation made use of CLSI EP6-P2.<sup>15</sup>

**Table 5: Piccolo Dynamic Ranges**

Analyte	Common Units	SI Units
<b>CHOL</b>	20–520 mg/dL	0.52–13.5 mmol/L
<b>HDL</b>	15–100 mg/dL	0.39–2.59 mmol/L
<b>TRIG</b>	20–500 mg/dL	0.23–5.65 mmol/L

If the analyte concentration is above the measuring range (dynamic range), but less than the system range, the printout will indicate a “>” sign at the upper limit and an asterisk after the number, e.g. CHOL >520\* mg/dL. If lower than the dynamic range, a “<” will be printed with an asterisk, e.g. CHOL <20\* mg/dL. For values that are grossly beyond the measurement range (system range), “~~~” will be printed instead of a result. Any time “~~~” appears on a printout, collect a new sample and rerun the test. If results for the second sample are suppressed again, please call Abaxis Technical Service.

**Sensitivity (Limits of Detection)**

The lower limit of the reportable (dynamic) range for each analyte is: cholesterol 20 mg/dL (0.52 mmol/L); HDL 15 mg/dL (0.39 mmol/L), and triglycerides 20 mg/dL (0.23 mmol/L).

**Precision**

Precision studies were conducted using CLSI EP5-A guidelines<sup>16</sup> with modifications based on CLSI EP18-P<sup>17</sup> for unit-use devices. Results for within-run and total precision were determined using two serum samples. The studies made use of data from four analyzers, two lots of discs, and two different sites across 10 days for both samples. The results are summarized below.

**Table 6: Precision**

Analyte	Sample Size	Within-Run	Total
<b>Total Cholesterol (CHOL)</b>			
<u>Serum 1</u>	N = 160		
Mean (mmol/L)		5.79	5.79
SD		0.08	0.15
%CV		1.3	2.6
<u>Serum 2</u>	N = 160		
Mean (mmol/L)		5.24	202.2
SD		0.08	0.11
%CV		1.5	2.2
<b>HDL</b>			
<u>Serum 1</u>	N = 160		
Mean (mmol/L)		1.43	1.43
SD		0.04	0.05
%CV		2.6	3.5
<u>Serum 2</u>	N = 160		
Mean (mmol/L)		0.98	0.98
SD		0.03	0.04
%CV		3.5	4.3
<b>Triglycerides (TRIG)</b>			
<u>Serum 1</u>	N = 160		
Mean (mmol/L)		2.34	2.34
SD		0.05	0.06
%CV		2.3	2.6
<u>Serum 2</u>	N = 160		
Mean (mmol/L)		1.85	1.85
SD		0.02	0.03
%CV		1.1	1.5

This data indicate that all three assays meet the NCEP precision criteria.<sup>2,3,4</sup>

### Correlation – Venous Sampling

Serum samples were collected and assayed by comparative methods. All test results were generated at a field site. The samples were chosen in order to provide a distribution of values using CLSI EP9-A2 guideline as a target for each analyte.<sup>18</sup> Representative correlation statistics are shown in Table 7.

**Table 7: Correlation with Comparative Methods\***

Assay	Correlation Coefficient (r)	Slope	Intercept	SEE	N	Sample Range (mmol/L)	Comparative Method
<b>Cholesterol (CHOL)</b>	0.997	1.079	-17.1	4.5	174	2.98-342	Bayer Cholesterol Assay on Hitachi 917
<b>HDL</b>	0.965	0.851	8.3	3.9	166	0.6-2.5	Roche HDL-C plus on Hitachi 917
<b>Triglycerides (TRIG)</b>	0.999	0.983	8.2	4.4	172	0.43-487	Bayer Triglyceride Assay on Hitachi 917

\*Data was generated using equivalence device Piccolo Blood Chemistry Analyzer

**Table 8: Calculated Recovery of the Abaxis Lipid Panel Assays – Venous Sampling**

Assay	Predicate Device Concentration mmol/L	Calculated Piccolo Recovery from Linear Regression Data Above mmol/L	Bias mmol/L	% Bias
<b>Cholesterol (CHOL)</b>	5.18	5.15	-0.03	-0.5
	6.22	6.27	0.05	0.8
<b>HDL</b>	1.04	1.09	0.05	5.0
	1.55	1.53	-0.02	-1.7
<b>Triglycerides (TRIG)</b>	1.7	1.76	0.06	4.0
	2.26	2.32	0.06	2.5

**Correlation – Capillary Sampling**

Heparinized capillary whole blood samples were collected and tested in singlicate on the Piccolo Xpress chemistry analyzer. Matched venous plasma samples from the same subjects were tested in duplicate by Roche test methods. The capillary samples were tested in three non-laboratory settings and the data combined. The samples were chosen in order to provide a distribution of values using CLSI EP9-A2 guideline as a target for each analyte.<sup>18</sup>

Representative correlation statistics are shown in Table 9.

**Table 9: Correlation of Piccolo Xpress chemistry analyzer with Comparative Methods – Capillary Sampling**

Assay	Correlation Coefficient (r)	Slope	Intercept	SEE	N	Sample Range (mmol/L)	Comparative Method
<b>Total Cholesterol (CHOL)</b>	0.995	0.97	1.2	5.6	639	0.5 - 10.7	Roche Total Cholesterol on Cobas 6000
<b>HDL</b>	0.981	0.99	-1.6	2.7	559	0.5 – 2.4	Roche HDL-Cholesterol Plus 3rd Generation on Cobas 6000
<b>Triglycerides (TRIG)</b>	0.996	0.96	4.1	7.9	588	0.4 – 5.6	Roche Triglycerides on Cobas 6000

**Table 10: Calculated Recovery of the Abaxis Lipid Panel Assays – Capillary Sampling**

Assay	Roche Test Concentration mmol/L	Calculated Piccolo Recovery from Linear Regression Data Above mmol/L	Bias mmol/L	% Bias
<b>Total Cholesterol (CHOL)</b>	5.18	5.02	-0.16	- 3.0
	6.22	6.03	-0.19	- 3.3
<b>HDL</b>	1.04	0.98	-0.06	-5.0
	1.55	1.5	-0.05	-3.3
<b>Triglycerides (TRIG)</b>	1.7	1.67	-0.03	- 1.3
	2.26	2.21	-0.04	- 2.0

## Accuracy – Cholesterol Reference Method Laboratory Network (CRMLN) Certification

Accuracy of the Piccolo assays for total cholesterol and HDL cholesterol was established by the requirements of the “HDL Cholesterol Method Evaluation Protocol for Manufacturers” of the US National Reference System for Cholesterol, CRMLN (Cholesterol Reference Method Laboratory Network), October 2018 using serum. A key part of the CRMLN protocol is linear regression analysis of the Piccolo assays versus the reference methods. The accuracy of the Total Cholesterol Assay compared to Abell-Kendall reference method is indicated by the correlation coefficient ( $R^2$ ) of 0.996, slope of 0.972, and intercept of 7.2 mg/dL (0.19 mmol/L). An among-run (n=10) CV for the Piccolo Total Cholesterol Assay was determined to be 0.8%.

For the Piccolo HDL Assay compared to the HDL reference method, precipitation followed by Abell-Kendall cholesterol assay, the correlation coefficient ( $R^2$ ) was 0.986, slope of 0.968, and intercept of 2.1 mg/dL (0.05 mmol/L). An among-run (n=20) CV for the Piccolo HDL Assay was determined to be 1.9%.

The Piccolo Total Cholesterol and HDL assays are not certified by the Cholesterol Reference Method Laboratory Network (CRMLN), which is coordinated by the Centers for Disease Control. The requirements of accuracy and precision have been previously established to meet the accuracy and precision requirements of CRMLN.

### Results of Untrained User Study

An “untrained user” study was conducted in which participants were given only the test instructions and asked to perform testing of 3 lipid panels with blinded randomized samples. The samples consisted of serum pools prepared at three levels for each of the three analytes, total cholesterol, HDL cholesterol, and triglycerides. The participants were not given any training on the use of the test. A total of 63 participants were enrolled from 3 sites, representing a diverse demographic (educational, age, gender, etc.) population.

Tables below present the summary of the performance for each analyte.

#### Total Cholesterol

	Level 1	Level 2	Level 3
N	63	63	63
Mean	3.73 mmol/L	5.14 mmol/L	6.35 mmol/L
%CV	2.9%	2.3%	1.3%
Observed Range	3.16 – 3.99	4.82 – 5.75	6.13 – 6.6
Percent of Results in the Range ± 11.1%*	98.4% (62/63) 95%CI: 91.5% to 100%	100% (63/63) 95% CI: 94.3% to 100%	100% (63/63) 95% CI: 94.3% to 100%

\*This percent is based on the premise that one cannot distinguish properly between normal and abnormal values when errors are greater than one-quarter of the normal range. The range of (140 mg/dL -220 mg/dL) was considered.

#### HDL Cholesterol

	Level 1	Level 2	Level 3
N	63	63	63
Mean	0.76 mmol/L	1.15 mmol/L	1.53 mmol/L
%CV	3.3%	3.2%	2.0%
Observed Range	0.73 – 0.83	1.09 – 1.24	1.48 – 1.61
Percent of Results in the Range ± 15.0%	100% (63/63) 95%CI: 94.3% to 100%	100% (63/63) 95% CI: 94.3% to 100%	100% (63/63) 95% CI: 94.3% to 100%

**Triglycerides**

	Level 1	Level 2	Level 3
N	63	63	63
Mean	0.94 mmol/L	1.73 mmol/L	2.32 mmol/L
%CV	3.0%	1.5%	0.9%
Observed Range	0.87 – 1.08	1.67 – 1.85	2.27 – 2.37
Percent of Results in the Range ± 15.0%	100% (63/63) 95% CI: 94.3% to 100%	100% (63/63) 95% CI: 94.3% to 100%	100% (63/63) 95% CI: 94.3% to 100%

### 13. Symbols



Use By



Catalog Number



Batch Code



In Vitro Diagnostic  
Medical Device



Consult Instructions  
for Use



Manufacturer



Do Not Reuse




X Number of Test  
Devices in Kit



Manufacturing  
Sequence



Serial Number

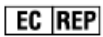

Authorized  
Representative  
in the European  
Community



Temperature  
Limitation



Caution

PN:  
Part Number



denotes conformity to specified  
European directives

## 14. Bibliography

1. Castelli, WP, et al. Lipids and risk of coronary heart disease. The Framingham Study. *Ann Epidemiol* 1992; 2:23-28.
2. National Cholesterol Education Program Expert Panel. Third report of the National Cholesterol Education Program (NCEP) Expert Panel on Detection, Evaluation, and Treatment of High Blood Cholesterol in Adults (ATP III). NIH Publication. Bethesda, MD: National Heart, Lung, and Blood Institute, 2002.
3. Executive summary of the third report of the National Cholesterol Education Program (NCEP) Expert Panel on Detection, Evaluation, and Treatment of High Blood Cholesterol in Adults (Adult Treatment Panel III). *JAMA* 2001; 285:2486-97.
4. Warnick, GR, et al. Impact of the third cholesterol report from the Adult Treatment Panel of the National Cholesterol Education Program on the Clinical Laboratory. *Clin Chem* 2002; 48:11-17.
5. Kayamori, Y, et al. Endpoint colorimetric method for assaying total cholesterol in serum with cholesterol dehydrogenase. *Clin Chem* 1999; 45:2158-2163.
6. Warnick G, Nauck M, Rifai N. Evolution of methods for measurement of high-density lipoprotein cholesterol: from ultracentrifugation to homogeneous assays. *Clin Chem* 2001; 47:1579-96.
7. Klotzsch SG, McNamara JR. Triglyceride measurements: a review of methods and interferences. *Clin Chem* 1990; 36:1605-1613.
8. Cole TG, Klotzsch SG, McNamara JR. Measurement of Triglyceride Concentration. In *Handbook of Lipoprotein Testing*, 2<sup>nd</sup> ed, Rifai N, Warnick GR, Dominiczak MH, eds. Washington, DC: AACC Press. 2000: 207-219.
9. Friedewald WT, Levy RI, Fredrickson DS. Estimation of the concentration of low-density lipoprotein cholesterol in plasma, without use of the preparative ultracentrifuge. *Clin Chem* 1972; 18:499-502.
10. Bachorik PS. Measurement of low-density-lipoprotein cholesterol. In *Handbook of Lipoprotein Testing*, 2<sup>nd</sup> ed, Rifai N, Warnick GR, Dominiczak MH, eds. Washington, DC: AACC Press. 2000: 245-263.
11. Clinical and Laboratory Standards Institute, CLSI (formerly, National Committee for Clinical Laboratory Standards, NCCLS). *Clinical laboratory waste management; approved guideline – second edition*. CLSI Document GP5-A2. Wayne, PA: CLSI, 2002.
12. CLSI GP41 7ED 2017: *Collection of Diagnostic Venous Blood Specimens*.
13. CLSI GP44-A4 2010: *Procedures for the Handling and Processing of Blood Specimens for Common Laboratory Tests; Approved Guideline—Fourth Edition*
14. Kroll MH, et al. Standardization of Lipoprotein Reporting. *Am J Clin Pathol* 2000; 114:696-702.
15. CLSI. *Evaluation of Linearity of Quantitative Analytical methods; proposed guideline – second edition*. CLSI Document EP6-P2. Wayne, PA: CLSI, 2001.
16. CLSI. *Evaluation of precision performance of clinical chemistry devices; approved guideline*. CLSI Document EP5-A. Wayne, PA: CLSI, 1999.
17. CLSI. *Quality management for unit-use testing; proposed guideline*. CLSI Document EP18-P. Wayne, PA: CLSI, 1999.
18. CLSI. *Method comparison and bias estimation using patient samples; approved guideline – second edition*. CLSI Document EP9-A2. Wayne, PA: CLSI, 2002.